

ARTRO

the ART of arthroscopy



ARTRO is a highly specialized electrocautery unit dedicated for bipolar and monopolar arthroscopic procedures.

ARTRO allows vaporization, cutting and bipolar coagulation in saline solution. It also offers monopolar cutting and coagulation modes which are necessary in some arthroscopic operations.

SDS system recognizes connected instruments and automatically selects the optimal working parameters. The operating modes and the output settings are selected automatically to the connected instrument.

Automatic regulation of output parameters guarantees fast and effective work.

CUT



VAPORIZATION



COAGULATION



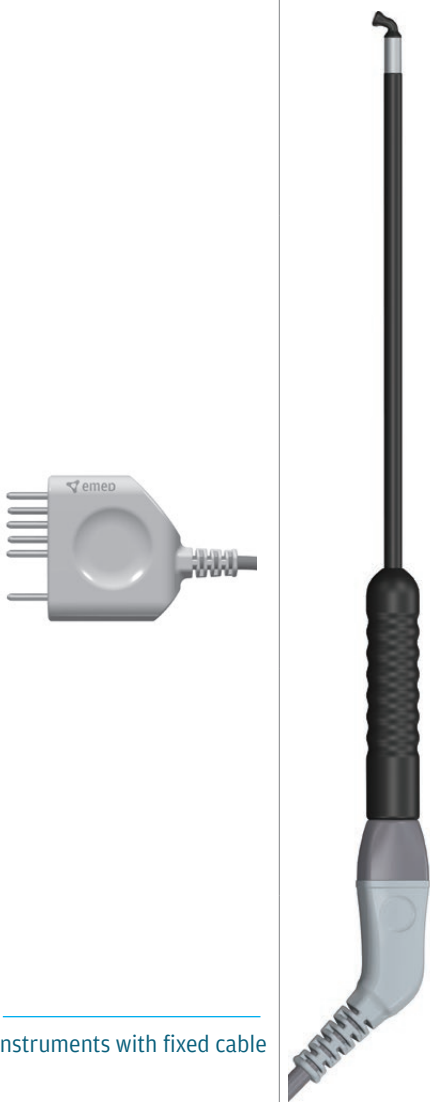






SmartDevice System

Advantages of the device

- versatile application range: monopolar and bipolar mode
- highly simplified operation
- color display with **touchscreen**
- advanced bipolar modes for use in saline solution
- automatic adjustment of the parameters - **SpectrumResult**
- scope protection
- recognition of the connected devices - **SmartDevice System**
- wide range of electrodes in different shapes and sizes
- reusable accessories to reduce the costs of the procedure
- dedicated **trolley** with convenient case for accessories



SDS bipolar arthroscopic electrodes

SDS plug	Bipolar electrode		REF		cable length	
 <p>instruments with fixed cable</p>			58S-010 58S-030	115 mm 170 mm	3 m	Bipolar needle electrode, angled 90°, 0.6 x 1.5 mm
			58S-020 58S-060	115 mm 170 mm	3 m	Bipolar needle electrode, angled 90°, 0.6 x 4 mm
			58S-040	115 mm	3 m	Bipolar button electrode, angled 90°, 1.5 x 4 mm
			58S-110 58S-140	115 mm 170 mm	3 m	Bipolar vaporization electrode, Phazer, convex, angled 70°, ball 2.4 mm
			58S-930 58S-920	115 mm 170 mm	3 m	Bipolar vaporization electrode, round, VAP

[100-610](#) Electrosurgical unit ARTRO



[100-302](#) Two-pedal footswitch, cable 5m, 6-pin plug



[080-061](#) TinyLine trolley for electrosurgical units



EMED SP. Z O.O. SP. K. | Ryżowa 69a | 05-816 Opacz-Kolonia | Poland
 Phone: +48 22 723 08 00 | email: export@emed.pl | www.emed.pl



Attention! This brochure does not replace the instructions for use! Consult instructions for use!

The information contained in this brochure is intended for health care professionals and medical equipment distributors.

Copyright © EMED All rights reserved. Any copying, distribution and publication in whole or in part without written permission of EMED is prohibited.

emed
 confidence in performance

REF 2023.01.K33.EN